PHONAK
life is on


A Sonova brand

Exchanging the earwax protection


Your hearing aid is equipped with an earwax protection to protect the speaker from earwax damage. Please check the earwax protector regularly and exchange it if it looks dirty or if the hearing aid volume or sound quality is reduced. The earwax protection should be replaced every four to eight weeks. To replace the earwax protector on the receiver, only use the wax filter tool that was recommended by your hearing care professional. Consult your hearing care professional about the availability of the wax filter tools.

Removing the earpiece from the speaker


1. Remove the earpiece from the speaker by holding the speaker in one hand and the earpiece in the other.

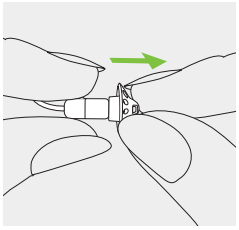2. Gently pull off the earpiece to remove.

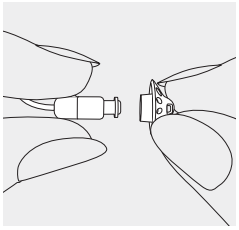
3. Clean the speaker with a lint-free cloth.

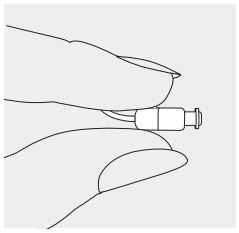
Exchanging the earwax protection with CeruStop

To replace the earwax protector on the receiver, only use the wax filter tool that was recommended by your hearing care professional.


1. Insert the removal side of the exchange tool into the used earwax protector. The shaft of the holder should touch the rim of the earwax protector.

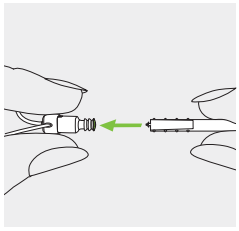
2. Carefully pull and remove the earwax protector straight out of the speaker. Do not twist the earwax protector when removing.

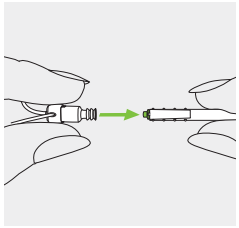3. To insert the new earwax protector, gently push the insertion side of the exchange tool straight into the hole of the speaker until the outer ring is perfectly aligned.

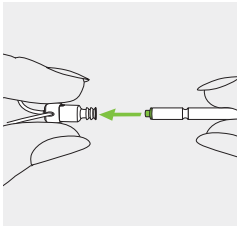4. Pull the tool straight out. The new earwax protector will remain in place. Now attach your earpiece to the speaker.

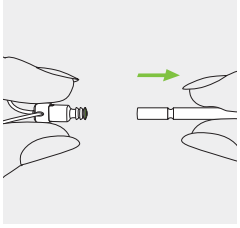Scan here for helpful hearing aid maintenance videos

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028-6848-03/V1.00/2024-08/hc © 2024 Sonova AG. All rights reserved.